

Caverject*impulse*[®]
Dual Chamber System
(alprostadil for injection)

A note on dosing:

- DO NOT use CAVERJECT IMPULSE more than 3 times a week
- Wait at least 24 hours between doses
- DO NOT change the dose without talking to your doctor
- DO NOT use any leftover CAVERJECT IMPULSE solution
- NEVER reuse the syringe or needle
- NEVER share or let anyone else use your CAVERJECT IMPULSE
- DO NOT use other medical treatments for ED while you are using CAVERJECT IMPULSE

Important Safety Information

CAVERJECT IMPULSE is used to treat erection problems in adult men. It is given by injection into the side of the penis. CAVERJECT IMPULSE should not be used by men with medical problems that can cause long-lasting erections—4 hours or more. Priapism requires immediate medical attention because permanent damage to the tissue of the penis and loss of potency may result. Tell your doctor about all of your medical conditions before using this drug. You must get medical attention right away if you still have an erection 4 hours after injection. This is dangerous.

The most common side effect of CAVERJECT IMPULSE is mild to moderate pain after injection. About one third of patients reported this pain but few stopped CAVERJECT IMPULSE for this reason. Men may get a small amount of bleeding at the injection site. Call your doctor if you notice any redness, lumps, swelling, tenderness, or curving of the erect penis. Lumps or curving of the erect penis including a condition called Peyronie's disease was reported in clinical trials with CAVERJECT IMPULSE. Please see enclosed product information.

At Home Injection Guide




**Once instructed by your doctor,
refer to this as needed.**

Caverject*impulse*[®]
Dual Chamber System
(alprostadil for injection)

Instructions to Prepare and Give an Injection of CAVERJECT IMPULSE

The more you use CAVERJECT IMPULSE, the easier it should get to use it. In time, you'll get into a routine. So be patient. Give yourself time to learn these steps. **REMEMBER:** if you have questions, ask your doctor.

The parts of your injection device:

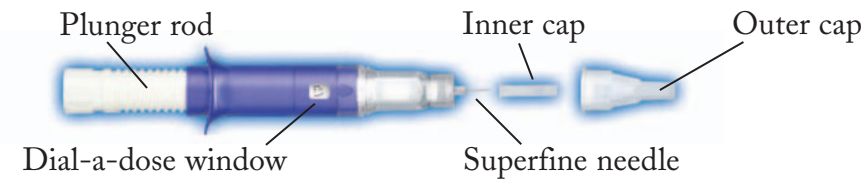

1

PREPARING FOR AN INJECTION



Wash your hands well with soap and water. Dry them with a clean towel. Take these items out from the tray and put them on a clean, flat, dry surface:

- Injection device
- Needle
- Alcohol swabs

As with any product, check the expiration date. Make sure that the date you inject the product is before the expiration date shown. Put the needle on the device:

- Use 1 of the alcohol swabs to clean the rubber tip of the syringe
- Peel the paper lid from the needle
- Press the needle onto the tip of the device. Turn it to the right until it is firmly in place
- Take the outer cap off the needle

2

MIXING THE SOLUTION



- Hold the device with the needle pointing up. The white plunger rod should be all the way out
- Turn the plunger rod slowly to the right until it stops. This will mix in the powder. Turn the device upside down 4 or 5 times to mix the solution evenly
- Carefully take the inner cap off the needle. Keep the device upright. Push the plunger rod in as far as it will go to get rid of the bubbles
- And pull the plunger back out again

Do not use if solution is cloudy or contains small pieces.

3

SELECTING THE DOSE



The window on the side of the device shows the dose:

- Hold the end of the plunger
- Turn it slowly to the right to dial-a-dose your doctor has prescribed. You may have to turn the rod a couple of times to start to see numbers
- Stop when you see your number. If the number you see is not the right number, keep turning the plunger rod slowly to the right. Set device down on a flat surface so the needle does not touch the surface

4

CHOOSING AND PREPARING AN INJECTION SITE



You should sit up straight or lie back slightly when you inject CAVERJECT IMPULSE. CAVERJECT IMPULSE is injected into the spongy tissue of the penis. It runs up and down your penis on the right and left sides.

- Choose an injection site on one side of the shaft of the penis
 - Make sure you avoid any veins you can see
 - Change the side of the penis you inject each time
 - Change the spot you inject each time
- Clean the injection site with a new alcohol swab
 - Do not throw this swab away. You will use it again
- Hold the head of your penis with your thumb and first finger
- Stretch your penis lengthwise along your thigh so the skin is tight

5

INJECTING CAVERJECT IMPULSE



Hold the penis firmly against your thigh to keep it from moving during the injection. Hold the device between your thumb and first finger.

- Push the needle into the injection site. Go through the skin and into the tissue as far as the needle will go until there is resistance
- Push the plunger rod in as far as it will go—until all of the medicine is injected
 - If the solution does not flow easily, move the needle a little bit and push again

When using a dose less than the full capacity, a small amount of liquid will remain in the device.

IMPORTANT: After you pull the needle out of your penis hold the alcohol swab on the injection site for about 5 minutes. Or until any bleeding stops.

- Carefully put the outer cap back on the needle
- Never share your medicine with anyone

6

CLEANING UP



The CAVERJECT IMPULSE syringe should never be reused.

- After injecting CAVERJECT IMPULSE, throw out the used syringe right away. Put it in a special container made for used syringes and needles. Like a hard plastic container with a screw-on cap. Or a metal container with a plastic lid labeled "Used Syringes." Do not use glass or clear plastic containers
- Learn the right way to throw out used vials, needles, and syringes. Be sure to ask your doctor or pharmacist what to do
- You can get these special containers from your doctor or pharmacist. Do not throw the container in the household trash or recycle bin. Always keep this container out of the reach of children